

MARS

BIO-MED PROCESSES INC



WE SAVE YOU MONEY!


Liberty **BOSS**

Amalgam Separator
ISO 11143 CERTIFIED 99.9%

1-866-594-3648

www.marsbiomed.com

PROPER INSTALLATION OF ANY AMALGAM SEPARATOR

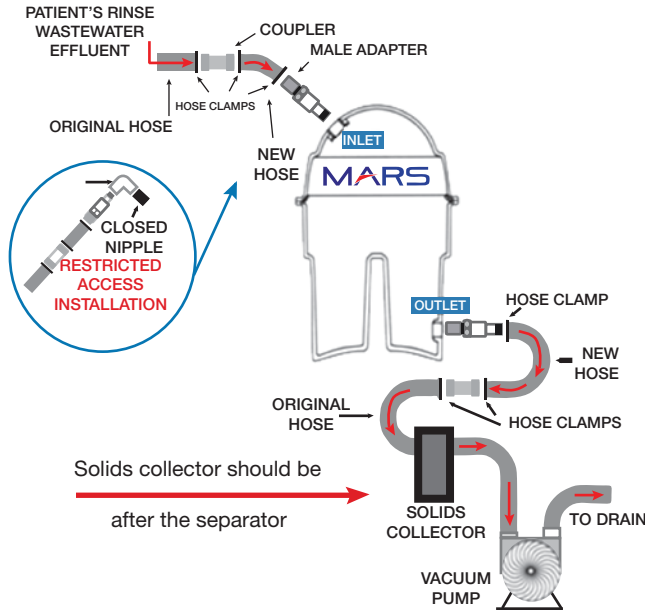

The **BOSS** will service 1-3 chairs for up to 3 years, **GUARANTEED.**

No hidden costs
No interim service
No messy filter changes

We encourage you to try and find a better guarantee, in writing.



Separator Service life

1 to 3 chairs: **Up to 3 yrs.**
 4 to 9 chairs: **Up to 2 yrs.**
 10+ chairs: **Up to 1 yrs.**

- LibertyBOSS separator uses an ion exchange media that traps mercury
- + 99% efficiency tested under actual operating conditions

How to Calculate Service Life

Although not an exact science, there is information required to give an educated estimate of service life;

1. How many total chairs
2. Number of Doctors
3. Number of Hygienists
4. Practice description (i.e. part time, average, busy)

With this information we will be able to give you an accurate estimate on your LibertyBOSS service life.

Performance Tested:

American Dental Association tested at 99.4% (2012)
 ISO 11143 Certified at 99.9%

PowerLINZ Evacuation Line Cleanser

All evacuation line cleaners that are used with amalgam separators should have a neutral pH. Highly alkaline or acidic cleaners will release mercury and decrease the efficiency of the separator. In time, non-neutral cleaners will also negatively impact your evacuation lines and your vacuum pump system.



PowerLINZ is fast and effective. Our formulation controls the rapid bacterial build-up inside the suction lines and the amalgam separator.

Everything you need in a line cleaner:

- Concentrated (1 bottle makes 125 liters)
- Cost-Effective
- Easy-to-Use
- Convenient Premeasured Pump
- Non-Corrosive
- Non-Foaming
- Antimicrobial
- Neutral pH
- Chlorine-free
- 100% Biodegradable

FOR MORE INFORMATION CALL:

1-866-594-3648

www.marsbiomed.com

