


PROPER INSTALLATION OF ANY AMALGAM SEPARATOR

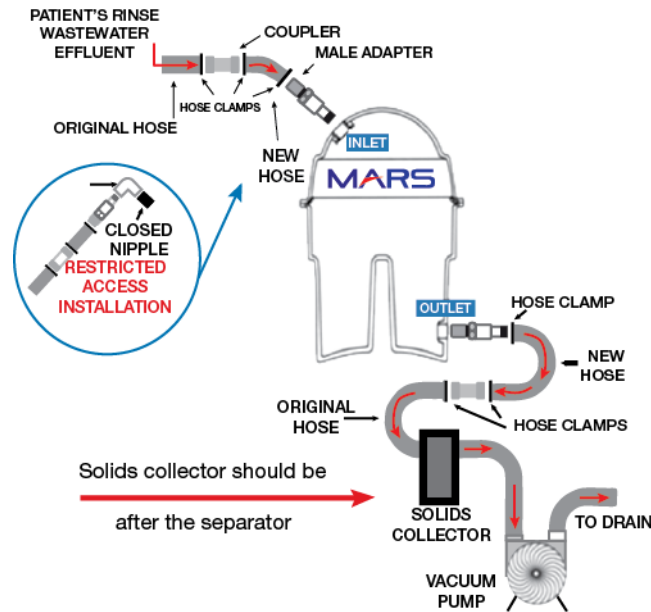

As seen in the above diagram, the solids collector (Evac-u-trap) is positioned **AFTER** the amalgam separator. This is because the discharge from the amalgam separator is safe for release into the sewer system. For your practice, this means that the material captured by the solids collector positioned after the amalgam separator can be disposed of as standard waste or flushed down the drain.

Placing the solids collector in front of the separator may appear to make some separator filters last longer, however, it also means that all waste collected by the Evac-u-trap should be considered toxic and must be disposed of in accordance with M.O.E. Toxic waste regulations.

M.A.R.S. Considers this unacceptable as it means the doctor is now responsible for two toxic waste disposals when only one is required.

We see far too many practices with incorrect separator installations that cost the practice thousands of dollars per year. To correct this problem and save money, call M.A.R.S. Bio-Med Processes Inc. at 905-723-9600 or 1-866-594-3648.